



# UROLOGY

REUSABLE INSTRUMENTS



+254 742 844 419






info@larioxhealth.com



www.larioxhealth.com

## WORKING ELEMENT MONOPOLAR & BIPOLAR (COMBO)


- Monopolar Resectoscope can be used in complex intrauterine operations without any difficulty
- Bipolar Resectoscope are particularly well suited for infertility patients
- Bipolar resection benefits from generation of plasma to achieve outstanding cutting power
- By using HF current several operations like hysteroscopic endometrium ablation, Myoma resection, adhesiolysis or septum dissection can be performed
- An electrolyte free irrigation fluid is used as distension media

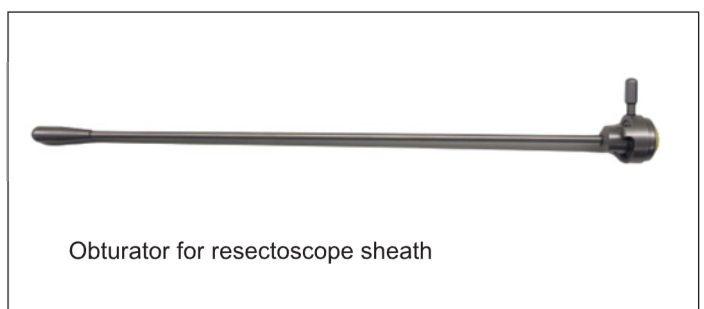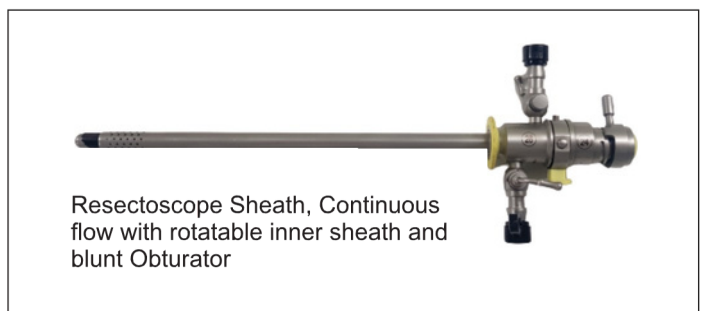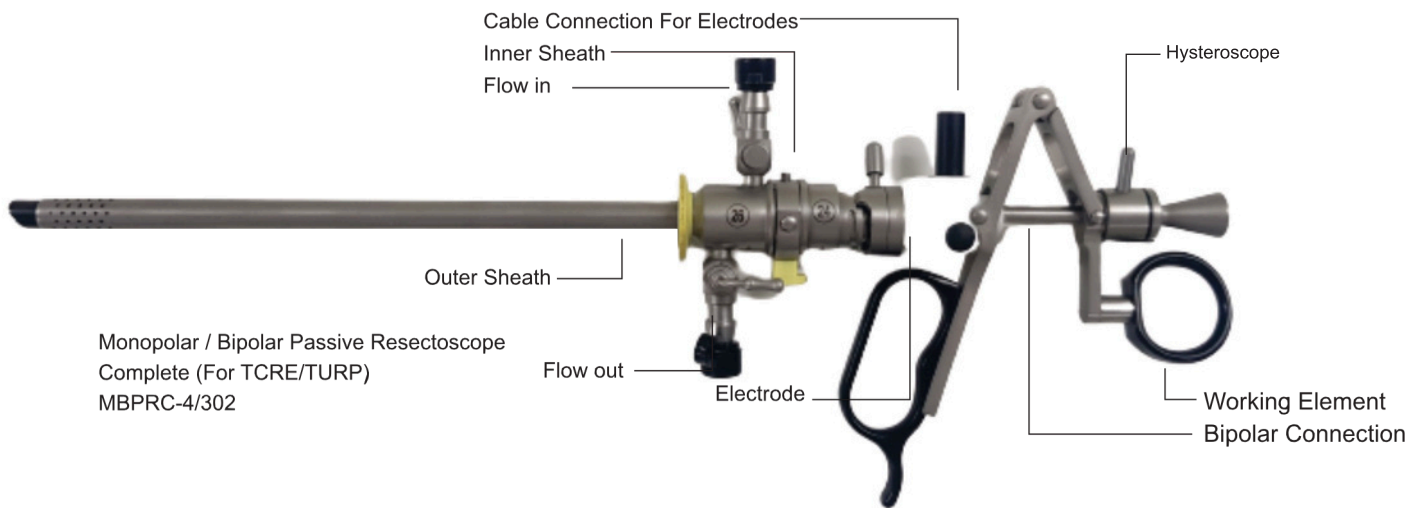

## STONE PUNCH







Bridge Single Channel  
and Double Channel

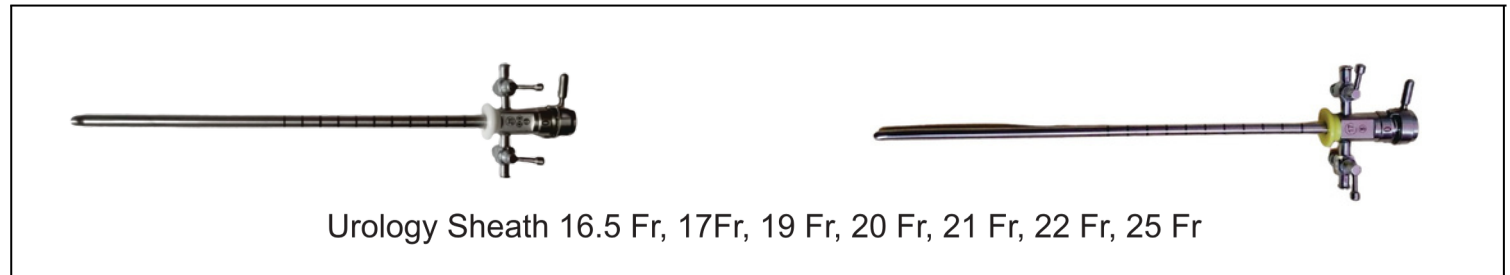

## ELECTRODES

HF- Electrodes, monopolar, single stem (for 24 fr)

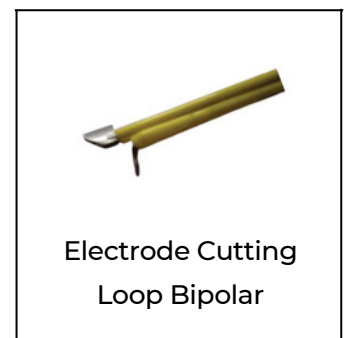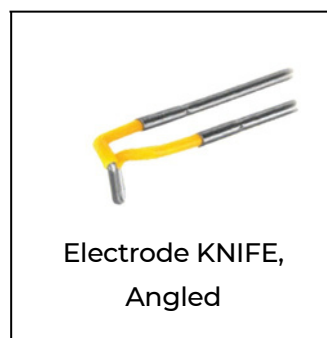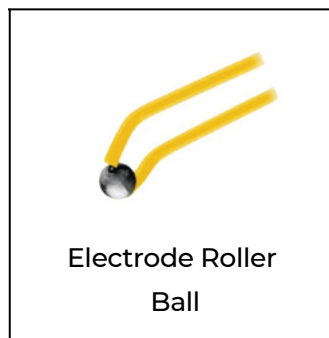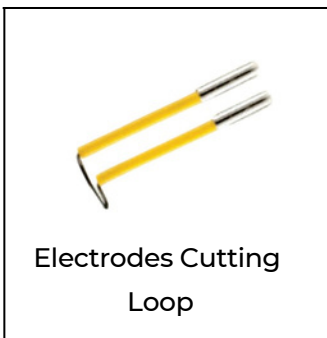



# CY series

## Flexible Video Cystoscope

### Suction function

One-button control to maintain clear view.  
 Unique anti-drop design to prevent the suction valve from falling off.



\*Patent NO.: ZL202022141662.4

### Multi-function buttons

- ① Zoom in/out
- ② Capture/record
- ③ Freeze/release



### Angle lock design

Lock or release the angle of the distal end for accurate diagnosis and treatment.



\*Patent NO.: ZL202022141661.X

### Tapered tip


Tapered tip design for easier insertion.




### Product Parameter

Model	Insertion Tube	Working Channel	Working Length	Bending (Up/Down)	Rotation (Left/Right)	Field of View	Depth of field	Suction function
CY-1356	15.6Fr	6.6Fr	380mm	220°/130°	120°	120°	3-100mm	Yes
CY-1355	15.6Fr	6.6Fr	380mm	220°/130°	120°	120°	3-100mm	No


### Image Processor




◀ HV-3101  
Image Processor




10.1"  
Touch screen




1280\*800  
HD Resolution




HDMI+ AV  
output



Detachable & rechargeable battery  
supporting 4 working-hour



Detachable  
SD card



Portable  
design



# UR series

## Flexible Video Ureterorenoscope

### Secondary bending

Bending U/D 275°  
Secondary bending to improve operation flexibility




### Tapered tip

Tapered tip design for easier insertion



### Rotation

Improve operation flexibility and accuracy

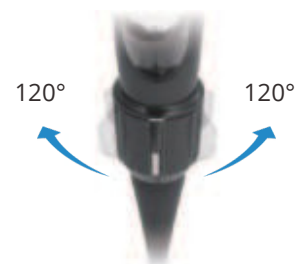

\*Patent NO.: ZL201621392021

### Multi-function buttons

- ① Zoom in/out
- ② Capture/record
- ③ Freeze/release



## Product Parameter

Model	Insertion Tube	Working Channel	Working Length	Bending (Up/Down)	Rotation (Left/Right)	Field of View	Depth of Field
UR-1328	7.5Fr	3.6Fr	670mm	275°/275°	120°	120°	3-50mm

## Image Processor



10.1" Touch screen



1280\*800 HD Resolution



HDMI+ AV output




Rechargeable battery supporting 4 working-hour




Detachable 64G SD card



Por table design

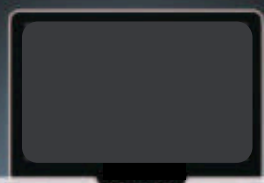

# Flexible Video Cystonephroscope

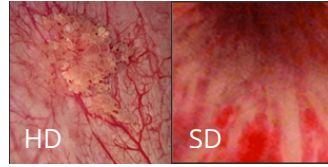

## HD Image

HD image gives more details of bladder.

## Control Button

- ① Zoom
- ② Capture/Recording
- ③ Freeze



## CMOS + LED Technology

Fiber-free design makes it more durable.



Insertion tube: 15.6Fr

Working channel: 6.6Fr

## Portable

Equipped with image processor HV-3101  
Support external monitor



## Application

Bladder examination, therapy and kidney calyx observation during PCNL.

Model	Insertion Tube	Channel	Working Length	Bending (Up/Down)	Field of View	Depth of Field
CY-1356	15.6Fr	6.6Fr	380mm	210°/130°	90°	3-50mm